



Exceptional stability and
durability for children on the go.



Durability and Stability

The Endolite Child's 4-bar Knee is designed for an active lifestyle. Featuring a polycentric design for stability during stance and an adjustable extension assist spring to provide swing control. It is a durable, lightweight and low maintenance knee designed for energetic children.



“The knee is very easy to work with. Adjustments are fast and easy. You can't go wrong with the tried and true 4-bar design.”

Scott, CPO

Versatility

Available with a choice of proximal adapters including a Lamination Socket Adapter Kit, allows for A-P shift or an M30 Threaded Adapter Kit, both are ideal where a minimal build height is required.

The knee features a 25mm diameter distal clamp, easy replacement of the pylon length, allows the system to grow with the child.

The 4-bar knee can be used with a variety of feet including the new Endolite Mini BladeXT. This multi-use children's foot offers outstanding stability, ground compliance and improved comfort through the split toe and traction sole, while its C-shaped toe spring is primed for optimal energy response.

Paired together the 4-bar Knee and Mini BladeXT make an unstoppable combination for children on the go.




Features

- Polycentric knee for stance stability
- Adjustable extension assist spring
- Proximal and distal attachment options
- Durable, lightweight design

Technical Information

- **Maximum User Weight:** 60kg (132 lbs.)
- **Activity Level:** 1-4
- **Component Weight:** 335g (12oz)
- **Build Height:** 105mm
- **Warranty:** 24 months

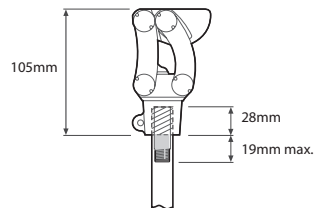

Part Number

Knee only - 239133
Knee w/pylon - 019245

Proximal Attachments

4 Prong Lamination
Anchor - 189626



M30 Threaded Adapter
Kit - 189726



Lamination Socket
Adapter kit - 189126



Distal Attachments

25mm Pyramid Pylon
Adapter - 409099



25mm Tube Clamp
Adapter - 189526




25mm Carbon Fiber
Pylon 180mm - 330129

25mm Carbon Fiber
Pylon 254mm - 330126



“I like the look of this one, its so small and simple looking. Feels very light weight, which is important for me and my short limb, looks cool like something in Star Wars.”

Braulo

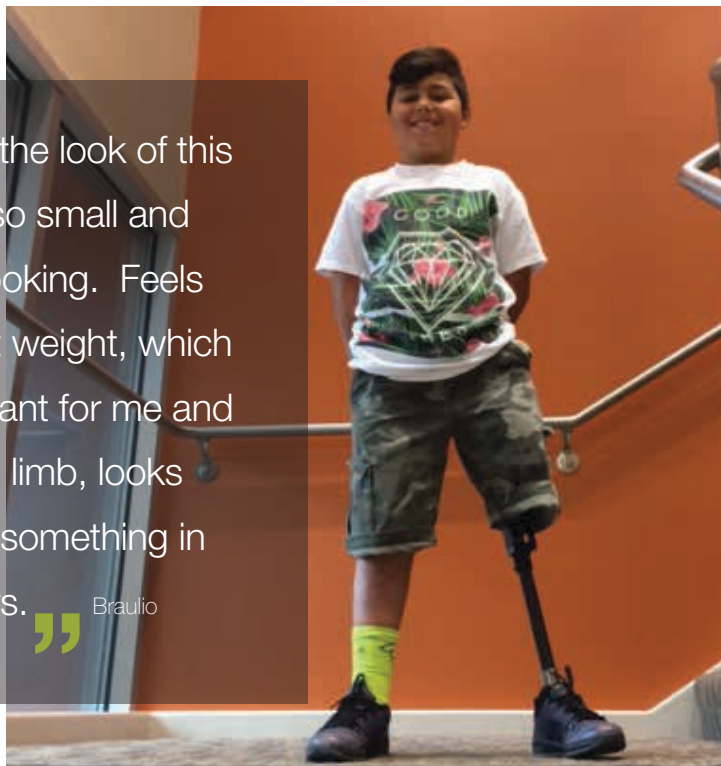

Follow Endolite:

 @EndoliteUS

 /EndoliteUS

endolite.com

info@endolite.com



Endolite is part of the Blatchford Group. Endolite Corporate Headquarters, 1031 Byers Road, Miamisburg, Ohio 45342
t: 800.548.3534 f: 800.929.3636

CE An annual visual inspection is recommended. Check for visual defects that may affect proper function. Maintenance must be carried out by competent personnel. Please check with your clinician before carrying out any new activities of daily living that you have not received training for with your Child's 4-Bar Knee.

Child4Bar US Iss1 01/17. Information correct at time of print.